

Best solution for
Your skin

BROW LIFTING

CHEEK LIFTING

WRINKLES IN THE
CORNER OF MOUTH

LIFTING NECK AND CHIN




Himedical!

#712 Medical Device Complex Center, 200 Gieopdosi-ro, Wonju-si,
Gangwon-do, Republic of Korea (26354)

Tel +82-33-747-6231 Fax +82-33-748-6231 E-mail sales@himedical.net
www.himedical.net

HIFU PLUS PREMIUM RF FACE LIFTING

- Non-invasive and Non-Surgical lifting
- Skin Tightening and Rejuvenation
- Eyebrow Lifting
- Nasolabial Fold Reduction

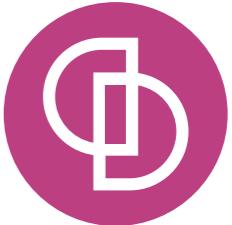

DUOLIF
DOUBLE CARE LIFTING

Double Lifting with HIFU
and PREMIUM 7MHz RF

DOUBLE CARE LIFTING

HIFU Mode


HIFU, High Intensity Focused Ultrasound is delivered into 3.0mm depth of skin fractionally and creates thermal coagulation without surface skin damages. The thermal energy delivered to SMAS layer helps skin tightening, pore improvement and collagen remodeling.

RF Mode

Multi-polar RF energy stimulates all layers of the skin, from epidermis to deep dermis and improves cell activation and blood circulation for skin lifting and tightening.

HIFU & RF Mode

Multi-polar RF energy is precisely focused at the coagulation spots that HIFU already delivered. RF energy remodels collagen at the targeted tissues for skin lifting and tightening and thermally stimulates them with more efficiently.

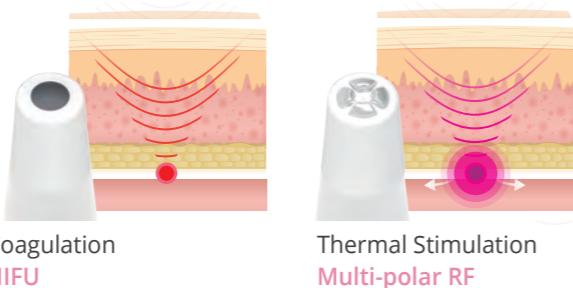

PREMIUM 7MHz DUOLIF

7MHz has low body impedance and enable to use low voltages so that it effectively generates deep thermal energy into the skin and reduces the risk of burn and electric shots compared to conventional RF devices in the market.

"Double Care Lifting! Using 10minutes a day!"

Effective


Improving the appearance of sagging skin around face, jowl and neck.

Proven results

It stimulates the deep structural support layers of the skin(3.0mm) without disturbing the surface of the skin, RF energy helps collagen regeneration for skin tightening in a short time and more efficiently.

Safety

HIFU and 7MHz RF technology delivers thermal energy into the deep inside without surface damage of the skin such as electric shots and burns.

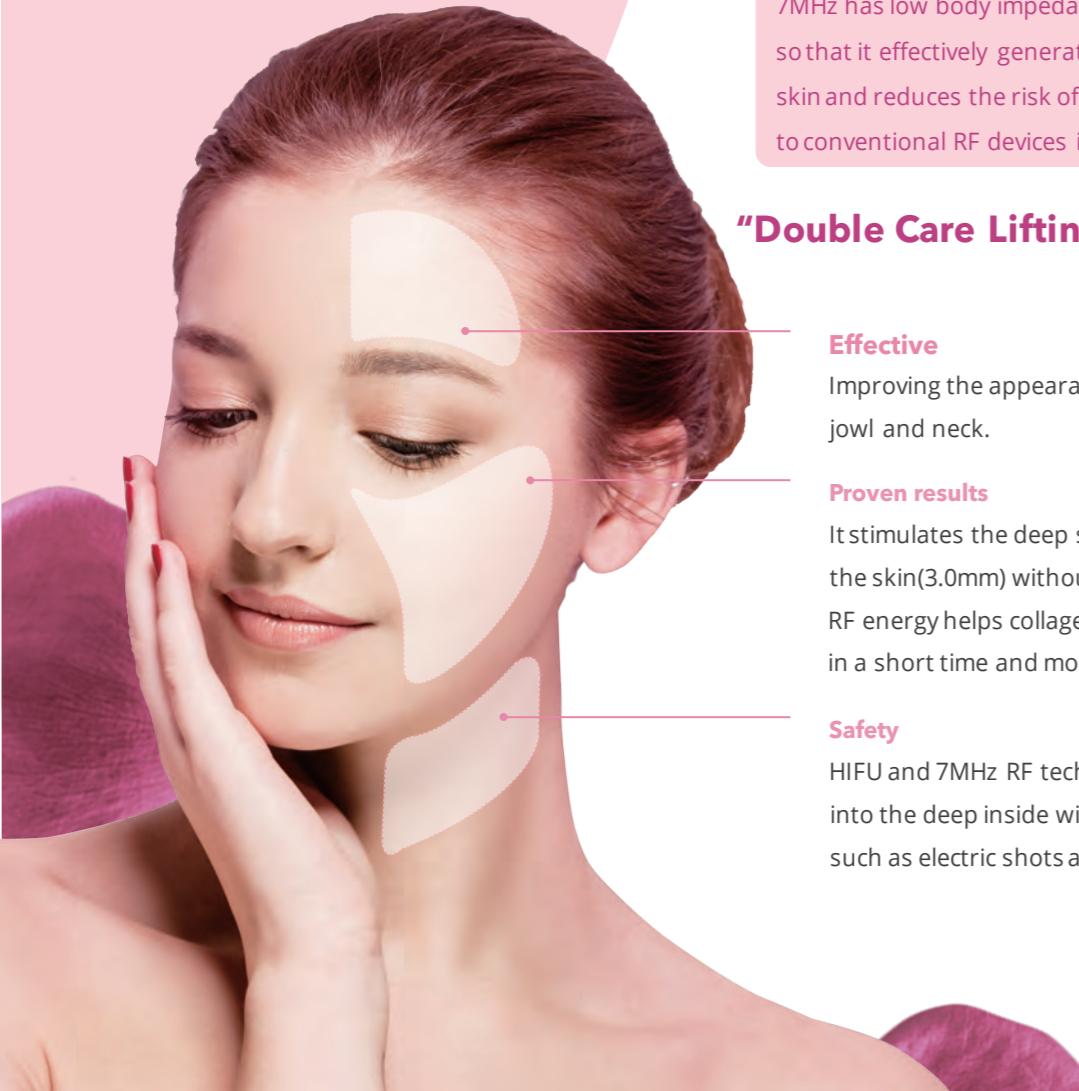

Bluetooth™



DESIGN AND SAFETY

Ergonomic Design

- Beautiful and compact size for easy using
- Simply changing cartridge with pressing the button
- Good Design with Hidden LED Display

Safety

- Safety and Efficiency with KC, CE certificate
- Avoiding double shots at the same spot with Gyro sensor(HIFU)
- RF is only emitted when the electrodes are contacted to the skin



ADVANTAGES

HIFU plus PREMIUM 7MHz RF

HIFU and RF mode are available in one main body. Transmitted HIFU penetrates to the depth of the SMAS layer without any damage to the epidermis.

The Multi-polar RF energy is specially focused on the coagulated spots that the HIFU energy emitted. The RF thermal energy activates the coagulation spots and activates the collagen remodeling and wound healing process more effectively.

Wireless Charging

Simply wireless charging when the device is placed into the charging cradle.

Cartridges

Using 60,000 shots for HIFU cartridge. It can be used about 6 years without exchanging the cartridge. Simply changing HIFU or RF cartridge by pressing one touch type button.

Hidden LED Display

Easy and intuitive recognition on LED display for mode, power level, cartridge information, remaining battery and charging level.

