



What is JILLSOF?

JILLSOF is an automatic drug delivery system optimized for treatment of vaginal dryness.



Development Background

There have been consistent needs and demands to improve inconvenient usability of syringe type HA that leads to long operation time during HA filler procedure in gynecology.

Dermashine(drug delivery system) and 32G-9 Pin multi needle, the representative products of Huons Medical, are the world's best-selling products which have been sold the most among competing products in dermatology. We are confident of our technology and clinical experiences.

We focused on resolving existing problem caused when doctors are using HA filled syringe, so JILLSOF is designed to provide doctors with simple and fast HA filler procedure.


Product Advantages

32 Gauge – 4 Pin Multi Needle (V-Needle®)

- ✓ The micro thin needle of 32 gauge relieves patient's pain
- ✓ 4 Pin multi needle reduces operation time sharply.

Auto - Sensing Suction

- ✓ The auto-sensing suction function not only helps injecting HA into skin safe and reduces loss of solution by making needles to sticks to skin surface completely.

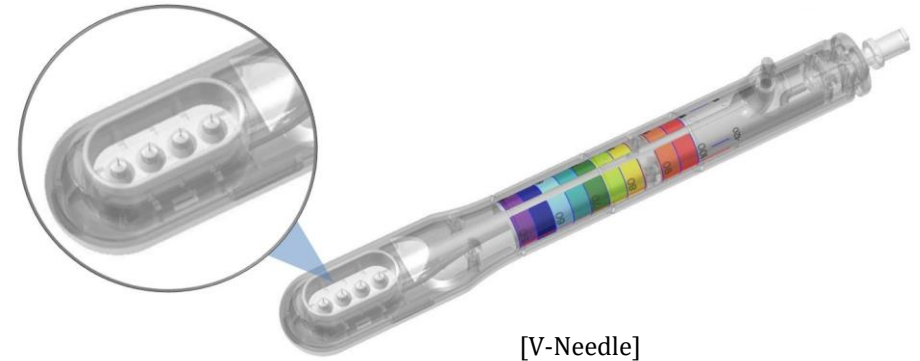

[V-Needle]



Frequent Asked Questions

What sets JILLSOF apart from other competitors?

- ✓ JILLSOF and multi-needle (V-Needle) are the rare products designed to inject HA into vagina safe and fast for treatment of vaginal dryness. There is no competitors capable of comparing with.

How long is operation time?

- ✓ It can vary in doctor's skill and procedure way. However, HA filler procedure using JILLSOF can shorten operation time by four(4) times than using general (pre-filled) syringe because the multi needle (V-Needle[®]) is composed of four(4) thin needles.

Is it painful?

- ✓ The multi needle(V-Needle[®]), which is composed of 32 gauge of micro-thin needles, reduces patient's pain.

How long does treatment effect last? (How long is operation interval)

- ✓ It is recommended to have a operation per 3 month ~ 6 month.



Product Summary

► Indication

- Vaginal dryness

► Advantages

- Rare item for treatment of VD
- 32 Gauge - 4 Pin multi needle
- Auto - sensing suction

Item	Drug Delivery System
Model Name	Panace-DS-52
Operation Mode	Dose/ Continuous-Normal/ Continuous-Slow
Infusion Rate	Min. 10 – Max. 140 Shot (per a syringe)
Syringe Option	1ml, 1.5ml, 2ml, 2.5ml, 3ml, 4ml, 5ml
Vacuum Pressure	Min. 200 mmHg – Max. 610 mmHg
Dimension	379mm x 328mm x 178mm
Weight	3.2 Kg (7lb)
Component	Main unit, injector, adapter, power code, foot switch,

